

SIRCHIE®

FINGER PRINT LABORATORIES

IODINE FUMING PRODUCTS

Fumette Iodine Fuming Guns

The **DF2016** consists of six (6) compact fuming guns. Its design permits fuming both horizontal and vertical surfaces. It is also excellent for fuming small objects with iodine vapor because it gives the user total control of the fuming operation. A hinged cap allows iodine containment when not in use.



No. DF2016 in use fuming document.

Iodine Print Enhancers



The **DCA16** consists of six (6) clean, convenient and easy-to-use enhancers. Each crushable glass ampoule contains 1ml of benzoflavone reagent housed in a plastic applicator tube. Unlimited storage life and completely disposable.

IO-FUME™ Instant Vaporizer

The **LPF1006** eliminates the need to create heat, developing latent prints quicker and easier! Use with any SIRCHIE fuming chamber, especially disposable chambers such as FR150 Insta-FumeTent™, and FR175 and FR175G Shake-N-Fume™ Disposable Chambers.



No. LPF1006 Components

Iodine Ampoules

The **AMP2066** consists of six (6) 3g glass ampoules designed for quick processing in the field or lab using plastic or polyethylene bags. Process fresh latents up to 48 hours.



No. AMP2066

Iodine Crystal Ampoules

The **A211C** consists of six (6) 3g glass ampoules may be used with any fuming apparatus. They are absolutely safe for storage in kits without vapor loss. Plus, they're very convenient and clean.



No. A211C

ElectroFume™ Iodine Fumer

The **FUM250** kit serves as a quick, efficient means of developing latents on porous surfaces while providing significant safety for the technician by permitting distance from iodine fumes. It features a high-volume, light-weight, battery-operated air pump and is supplied with both Fumette Disposable Iodine Fuming Guns and Iodine Print Enhancers.

When fuming with FUM250, you must warm the iodine print enhancer by holding it in the palm of your hand.



No. FUM250 Kit and Components

CAT. NO.	DESCRIPTION	PRICE
FUM250	ElectroFume™ Iodine Fumer	\$99.69
DF2016	Fumette Disposable Iodine Fuming Guns, 6 ea.	\$28.75
DCA16	Iodine Print Enhancers, 6 ea.	\$23.25
LPF1006	IO-FUME™ Instant Vaporizer, 6 pk.	\$27.95
AMP2066	Iodine Ampoules, 6 ea.	\$32.50
A211C	Iodine Crystal Ampoules, 6 ea.	\$14.25

LATENT PRINT DEVELOPMENT

NINHYDRIN PRODUCTS

Ninhydrin Spray

No. 201C

Most widely used formulation of Ninhydrin in the world. It has developed prints over 2 years old in actual cases and 50-year-old prints in experimental cases.




Special Formula Ninhydrin Spray


No. NSI609

Ideal for developing latent prints on paper because it virtually eliminates ink from running.

Caution:
Flammable.

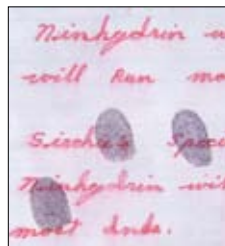

Ink ran on the document above when developed with acetone-based solvent, but using No. NSI609 Special Formula Ninhydrin on the document below to develop prints, the ink did not run.

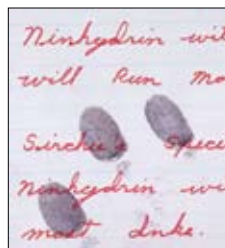


Cash register receipt developed visible prints when tested with NHT609 (close-up shown above).

Ninhydrin Spray with Acetone

No. 201ACE

Custom formulation of Ninhydrin in an acetone solution. It is used and recommended by experts in the latent print processing field for increased sensitivity.

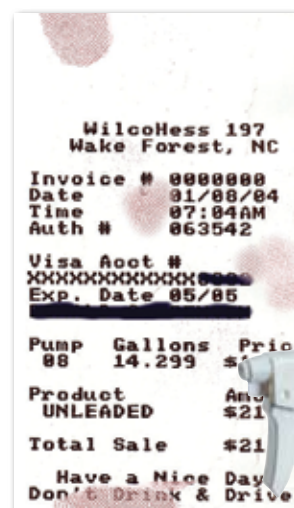

NHT609 developed latents on a gas pump receipt.

NinhydrinHT Spray

No. NHT609

Specifically formulated for recovering latents on carbonless and heat-transfer (thermal) papers, but also works on standard paper. Completely free of staining reagents, will not cause inks to run, and is non-flammable.



Ninhydrin Crystals

No. NRP01A

Ninhydrin Crystals are offered for those agencies that prefer to mix their own unique formulation.



Ninhydrin Fixative

No. NFS200

Prevent fading and preserve Ninhydrin developed latent prints after they've been photographed. In some cases, this fixative enhances and tints developed prints as well.



CAT. NO.	DESCRIPTION	PRICE
201C	Ninhydrin Spray, 16 oz. (473ml)	\$18.75
201ACE	Ninhydrin Spray with Acetone, 16 oz. (473ml)	\$17.50
NSI609	Special Formula Ninhydrin Pump Spray, 8 oz. (240ml)	\$31.95
NHT609	NinhydrinHT Pump Spray, 8 oz. (240ml)	\$49.00
NFS200	Ninhydrin Fixative Pump Spray, 8 oz. (240ml)	\$12.25
NRP01A	Ninhydrin Crystals, 25g	\$24.95